


# Jandel Universal Probe Guidance Unit



Needle  
Guidance  
Unit

The Universal Probe Needle Guidance Unit is supplied with a set of probe tips installed. Changing the probe tip spacing requires replacement of the Guidance Unit. Tip replacement can be performed by replacing the Guidance Unit, or just the needles can be replaced. Replacing the whole guidance unit is more complicated and unnecessary unless there is some damage to the current guidance unit. There are instructions for replacing the needles in the manual. Either method requires the ligaments to be unsoldered and then new ones soldered.

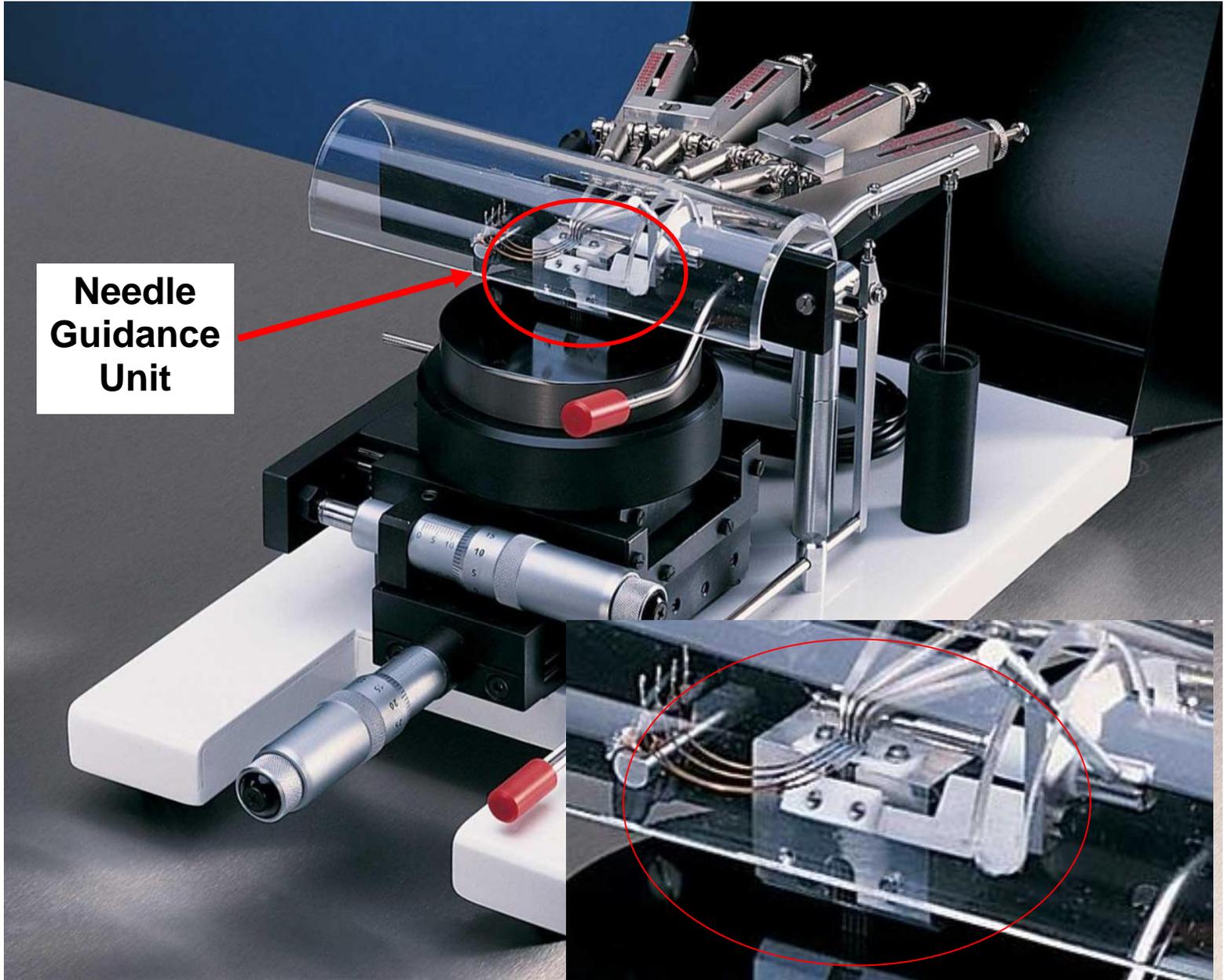